Saturn[®]

ST-BC1304

MASSAGER

APARAT DE MASAJ

GB

MASSAGER

Dear Buyer!

We congratulate you on having bought the device under trade name "Saturn". We are sure that our devices will become faithful and reliable assistance in your housekeeping.

Avoid extreme temperature changes. Rapid temperature change (e.g. when the unit is moved from freezing temperature to a warm room) may cause condensation inside the unit and a malfunction when it is switched on. In this case leave the unit at room temperature for at least 1.5 hours before switching it on.

If the unit has been in transit, leave it indoors for at least 1.5 hours before starting operation.

SAVE THESE INSTRUCTIONS

Turn massager on with the convenient three position switch that lets you choose heat only, massage only, or heat with massage.

Magnetize Infrared Heat with Massage (HM)

Infrared Heat with massage is most effective without the use of any attachments. Simply switch to "HM"(Infrared Heat and Massage) and enjoy the unique combination of penetrating infrared heat with a vibrating massage.

Magnetize Infrared Heat Only (H)

BEFORE using unit on Infrared HEAT, check the temperature of the infrared plate. You can control the amount of infrared heat you desire by switching unit from Off to On periodically.

Massage (M)

Four convenient attachments make it possible to choose the right massage for

each part of your body. Simply place the attachment you have selected over the heat plate of the massager. Use the massager wherever you need or want it-on arm, legs, back, thighs, neck, shoulders. Although the massager may be applied through clothing, it is most effective when used directly on the body.

Adjustable handle

You can adjust the angle of the massager for best massaging position. Simply press the circular button in the middle of the handle, hold it while you are adjusting the angle of the handle. Release the button upon reaching the desirable angle. You can adjust the angle up to maximum 90.

Roller Massager Head

Refreshing rollers with little stimunodes to soothe tired body.

Suction Head

Generates a slight vacuum, which has very beneficial tightening effect on your skin.

Body Massager Head

Specially designed double knodes direct soothing waves to shoulders, arm, wrist & legs.


Scalp Massager Head

Begins at the base of the scalp slowly moves up and over the crown of your head, then reverse direction for a invigorating massage.

Joint Massager Head

Move upwards with a firm, rotating motion relax all over body.

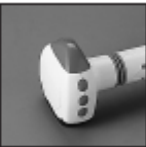Concentrator Head

Extended concentrating knodes provides gentle vibrating waves for wider body area such as legs, hip & back.

5. Do not reach for an appliance that has fallen into water, unplug immediately.

To reduce the risk of burns, fire, electric shock, or injury to persons:

1. Before using the appliance, make sure its voltage corresponds to the local power supply.
2. An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
3. Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
4. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
5. Keep the cord away from heated surfaces.
6. Never block the air openings of the appliance or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, etc.
7. Never use while sleeping.
8. Never drop or insert any object into any opening.
9. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
10. Do not use on or near eyes or other highly sensitive areas.
11. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
12. Do not carry this appliance by supply cord or use cord as a handle .
13. Do not use the vibrating massager directly on one concentrated area of the body. Vibrating massagers are meant to be used in a continuous motion, constantly changing the massage area.

14. To disconnect, turn all controls to the off position; then remove plug from outlet.

15. **CAUTION** - All servicing of this massager should be performed by authorized service personnel.
16. Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer. This appliance is designed for household use only.
17. Do not use outdoors.
18. Use heated surfaces carefully, may cause serious burns. Do not use over sensitive skin area or in presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.

1. Do not place in or drop massager or attachments in water or any liquid. Do not use while bathing or showering.
2. Always unplug massager immediately after use.
3. This massager should not be used by persons under the care of a physician without first consulting the physician. This device should not be used over swollen or inflamed areas or skin eruptions.
4. DO NOT continuously use for more than 20 minutes. Recommend to stop 20 minutes then restart to ensure LONGER LIFE of unit.

DIRECTIONS FOR USE

1. The massager is designed to work on ordinary household current only. It is not intended for commercial use.
 2. Unplug before cleaning. Clean only with a damp cloth. No other maintenance is required.
 3. Periodically straighten cord if it becomes twisted.
- Magnetic Infrared Lens
Built-in 3 pieces of magnets: The magnetic field is suggest that it helps in improving the blood circulation; thus enhance the result of massage.
Working of the infrared light

This light penetrates deep into the skin. The skin is therefore heated, so that the blood vessels are dilated and the flow of blood is stimulated. The blood brings with it the substances for rebuilding and nourishing the body tissues more rapidly whereas decomposition products are more quickly removed. The pores in the skin will also be further opened under the influence of the heat. These properties make it extremely suitable for therapeutic as well as for cosmetic applications.
Therapeutic applications
Infrared light can relieve the pain with rheumatic affections, muscular pain, we recommend that you consult your doctor first.
Cosmetic applications
Creams or other skin care agents are better absorbed under the influence of infrared light, and are therefore most effective.
The infrared heat will also help to keep your skin healthy and clean, because of its stimulating effect on the blood circulation.

Specifications

Power supply: 13W
Rated voltage: 220-230V
Rated frequency: 50Hz
Rated current: 0,06A

Set

MASSAGER.....1
 MESSAGE ACCESSORIES.6
 INSTRUCTION MANUAL WITH WARRANTY BOOK.....1
 PACKAGE.....1

ENVIRONMENT FRIENDLY DISPOSAL

You can help protect the environment!
Please remember to respect the local regulations: hand in the

IMPORTANT SAFEGUARDS

When using an electrical appliance, basic precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING

- To reduce the risk of electric shock:
- 1. Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
- 2. Do not use while bathing or in a shower.
- 3. Do not place or store appliance where it can fall or be pulled into a tub or sink.
- 4. Do not place in or drop into water or other liquid.

non-working electrical equipments to an appropriate waste disposal centre.

The manufacturer reserves the right to change the specification and design of goods.

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GB

INTERNATIONAL MANUFACTURER'S WARRANTY

The Warranty is provided for the period of 2 years or more if the Law on User Right of the country where the appliance was bought provides for a greater minimum warranty period.

All terms of warranty comply with to the Law on Protection of the Consumer Rights and are governed by the legislation of the country where the product was bought.

The warranty and free of charge repair are provided in any country where the product is delivered to by Saturn Home Appliances or authorized representative thereof, and where no restriction of import or other legal provisions interfere with rendering of warranty services and free of charge repair.

Cases uncovered by the warranty and free of charge repair:

1. Guarantee coupon is filled in improperly.
2. The product has become inoperable because of the Buyer's non-observance of the service regulations indicated in the instruction.
3. The product was used for professional, commercial or industrial purposes (except for the models, specially intended for this purpose that is indicated in the instruction).
4. The product has external mechanical damage or damage caused by penetration of the liquid inside, dust, insects and other foreign objects.
5. The product has damages caused by non-observance of the rules of power supply from batteries, mains or accumulators.
6. The product was exposed to opening, repair or modification of design by the persons, unauthorized to repair; independent cleaning of internal mechanisms etc. was made.
7. The product has natural wear of parts with limited service life, expendables etc.
8. The product has deposits of scale inside or outside of fire-bar elements, irrespective of used water quality.
9. The product has damages caused by effect of high (low) temperatures or fire on non-heat resistant parts of the product.
10. The product has damages of accessories and nozzles which are included into the complete set of product shipment (filters, grids, bags, flasks, cups, covers, knives, sealing rings, graters, disks, plates, tubes, hoses, brushes, and also power cords, headphone cords etc.).

On the warranty service and repair, please contact the authorized service centers of Saturn Home Appliances.

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GB

WARRANTY COUPON

**When purchasing the product,
please require its checking before you;
BE SURE
that the goods sold to you, are functional and complete
and that the warranty coupon is filled in correctly.**

This warranty coupon confirms the absence of any defects in the product you purchased and provides for free of charge repair of the product failed through the fault of the manufacturer throughout the period of warranty service and free of charge repair.

Failing the presentation of this coupon, in case of its improper filling in, infringement of factory seals (if any), and also in cases indicated in the warranty obligations, the claims are not accepted, and no warranty service and free of charge repair is made! The warranty coupon is valid only in the original copy with the stamp of trading organization, signature of the seller, date of sale and signature of the buyer.

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